

Fig 1

Features:

- 50Hz or 60Hz
- Medical grade stainless steel electrode contacts
- Protection against dust and water spray, compliant with IP64
- Dimensions 26.8mm x 18mm x 9.7mm (see Fig 1)
- Gain selection, 1-7 (see Fig 1)
- Gain and threshold can be changed using **biosim** software

Recommendations:

- Socket material should not be manufactured from carbon material as this can interfere with electrode performance
- For optimum performance ensure the ribbon cable is cut perpendicularly. The cable must be fully inserted to the block prior to assembly (see Fig 3)
- For optimum performance set gain to 4

Additional Kit Contents:

- Lamination and thermoplastic dummies (see Fig 2)
- Electrode cable (choice of 300mm or 600mm) and insulated displacement connector (see Fig 3)

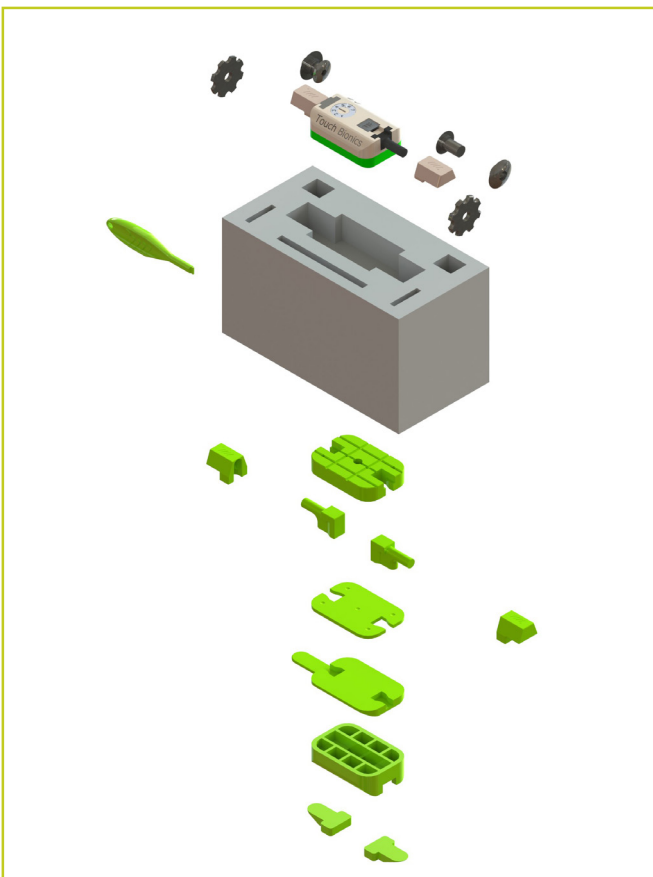


Fig 2

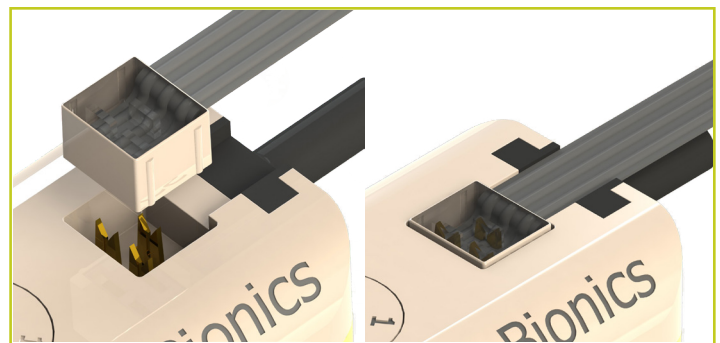


Fig 3

Part Number	Part Name
PL091050	Compact Electrode Kit 50Hz (300mm)
PL091060	Compact Electrode Kit 60Hz (300mm)
PL091127	Compact Electrode Kit 50Hz (600mm)
PL091128	Compact Electrode Kit 60Hz (600mm)