

Motion Control ETD



Heavy Duty Activity Package

Add the Motion Control titanium ETD with water resistant housing and sleeve to any i-limb purchase for heavy duty activities.


- Wide opening and rugged durability, combined with classic hook-type fingers
- High grip force of up to 25 lbs with proportional control allows users to effortlessly vary the grip force anywhere within the range
- Brushless DC motor technology is substantially more durable and rugged than existing motors in use
- Dependable in all conditions and environments due to the water resistant feature
- Titanium Hook Fingers offer close to twice the strength of typical aluminum fingers and are resistant to bending and scratching
- Interchangeable with i-limb hands utilizing quick wrist disconnect feature for heavy duty activities



Also available with **flexion wrist**

Features:

- Freedom and flexibility of natural wrist movement
- Smooth & natural motion which provides control in completion of daily tasks
- Enables shock and torque absorption
- Lock and unlock function via simple magnetic lock
- Available in black

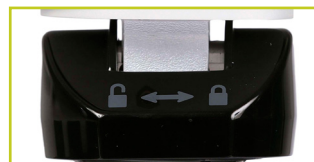

- Lock **flexion wrist** to position i-limb at +40°, +20°, 0°, -20°, -40°
- Unlock **flexion wrist** for springloaded flexibility range from +40° to -40°